



# Riser Control System

## Safety, functionality and ease of use...

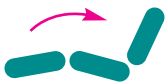

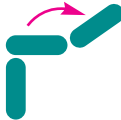



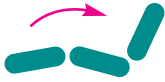


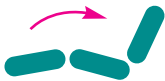


The Configura® Riser Control System provides a unique combination of safety and functionality for users with reduced cognitive ability.

The Configura® Riser Control System allows the therapist to disconnect selected chair functions, for example the riser, tilt or backrest recline functions, via a simple switch at the back of the chair.

For example, it may be beneficial for the user to be able to operate the tilt and recline functions, but unsafe to operate the riser function. The Riser Control System allows the riser function to be disabled on-site within moments via a sample switch.

The Configura® Riser Control System technology can be added to any chair at point of order.

## Safe, functional, simple...

<b>Setting 1 - Full function</b>	<b>1 Tilt-in-space</b> 	<b>2 Rise</b> 	<b>3 Recline (dual motor only)</b> 
<b>Setting 2 - no tilt-in-space</b>	<del><b>1 Tilt-in-space</b> </del>	<b>2 Rise</b> 	<b>3 Recline (dual motor only)</b> 
<b>Setting 3 - no rise</b>	<b>1 Tilt-in-space</b> 	<del><b>2 Rise</b> </del>	<b>3 Recline (dual motor only)</b> 
<b>Setting 4 - no recline (dual motor only)</b>	<b>1 Tilt-in-space</b> 	<b>2 Rise</b> 	<del><b>3 Recline (dual motor only)</b> </del>