

Invacare® | **Reliant 350 Stand Assist**



Fresh new look and
design improvements
for ease of use



For safe and comfortable assisted transfers

The *Invacare Reliant 350 Stand Assist* is a compact and maneuverable mobile lifter that provides safe and comfortable assisted transfers for those with limited mobility or rehabilitation needs. The *Invacare Reliant 350* can be used to transfer a client that cannot transfer unaided, but has some level of weight bearing ability, from one seated position to another. The robust design with a maximum weight limit of 158 kg (25 Stone) allows for additional support whilst transferring, preventing unnecessary strain on both patient and carer.



Yes, you can.®




Adaptable and Easy to Use

Reliant 350 features a multi-adjustable lift boom and a selection of sling mounting points to provide an adjustable and ergonomic lift angle. Large rear and small front castors help maneuver the **Reliant 350** into small or confined spaces and the base has a lever adjustment mechanism for added stability. In addition, it has a generous non slip footplate providing a stable surface for the feet, as well as a height-adjustable padded knee guard.

Safe and Reliable


The **Reliant 350** provides a range of safety features including an emergency stop button, emergency mechanical stop system, as well as an anti-entrapment device which stops the boom from lowering if it meets resistance. The battery recharge indicator is visual and audible. For transport and storage, the **Reliant 350** can be easily dismantled.

Features and Options

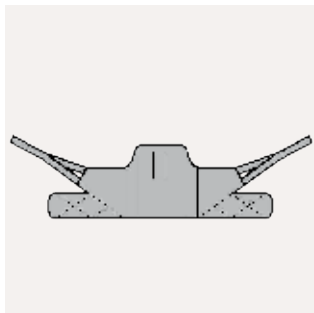

Stand Assist Sling

Designed to be used with a **Stand Assist** lifter, this sling features the latest innovative fabric for improved comfort and handling. Available in sizes; S, M, L & XL.

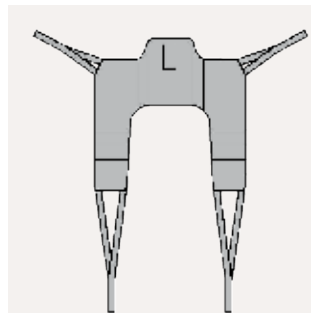

Transfer Stand Assist Sling

Designed for **Stand Assist** lifters, for those who require additional leg support, this sling provides a supportive and comfortable transfer. Available in sizes; S, M, L & XL.







Independent Charger

An off-board charger and spare battery allow for 24 hour use.

Technical data

Invacare Reliant 350					
	1670 mm	1020 mm	900 mm	940 mm	660 mm

Invacare Reliant 350				
	110 mm	100-200 lifts per full charge	49 kg	Max. 158 kg

INVACARE® Limited
Pencoed Technology Park - Pencoed -
Bridgend CF35 5HZ
Tel: +44 1 656 776 222
Fax: +44 1656 776 220
E-mail: uk@invacare.com
www.invacare.co.uk
Sales Order E-mail:
ordersuk@invacare.com

INVACARE® Ireland Ltd
Unit 5 - Seatown Business Campus
Seatown Road - SWORDS - County
Dublin
Tel. +353 1 8107084
Fax +353 1 8107085
Email: Ireland@invacare.com
www.invacare.ie

© 2011 Invacare® International Sarl
All rights reserved.

All information quoted is believed to be correct at time of print.
Invacare® reserves the right to alter product specifications without prior consultation.

Reliant 350 - UK - 03/2011 - Q99632

Frame colour

Product Code: RPS350-1E



Beige



Yes, you can.®