

Invacare® | **Robin**®



3 year warranty



Moving forward with dignity and care

The *Invacare Robin* ceiling hoist provides a unique and innovative method of transfer. The two strap design has no obstructive spreader bar; offering a safe and dignified transfer with optimal client and carer contact.



Yes, you can.®



Comfort

The two strap design gives unsurpassed stability during the transfer. The smooth soft start function makes positioning more gentle, and the two speed motor reduces the amount of time required for the transfer. The strap only configuration eliminates the need for a spreader bar of any size, adapting automatically to the weight and shape of the client.

Transfer with dignity


Personal contact is important. The absence of a spreader bar provides the client with a sense of freedom reducing the fear of being close to obtrusive parts of the lifting equipment. This enhances the transfer experience by introducing a unique dignified method of moving people.

High safety and easy to work with

The *Invacare Robin* has been conceived with safety in mind. In the event of an emergency situation, lowering, raising and stop functions are standard features built into the motor unit. The lifting mechanism offers an overload security

system, self adjusting strap inlet and centrifugal braking. Our patented strap system enables the carer to attach the sling to the hoist with the minimum of effort and the added security of a simple locking mechanism.


Features and Options

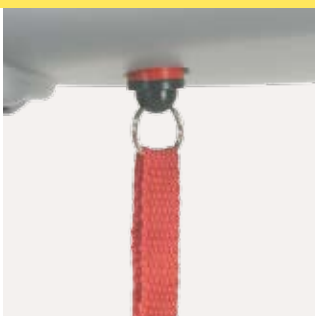

Emergency stop and lowering



Strap inlet self adjusts automatically to size of client

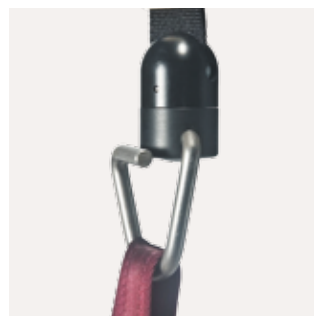

Hook self locks when weighted



Cable disengages to avoid damage if unit is pulled along via hand control

Features and Options



Robin
The unique strap system without the use of a spreader bar.



Robin Mover
With an integrated traverse motor.



Blue light illuminates when motor is in use, and a signal light flashes when battery is low



Charge motor through hand control when not in use



External emergency battery



Control- buttons on motor unit
For emergency use.



A scale which is easy to mount and use



Flexiscope
If a specific situation requires a spreader bar.










Stretcher, that offers unique stability without a leveller




① **Hand control**
Identified by symbols.
② **Robin Mover hand control**
With symbols for up/down and direction.

Technical data

					
Robin	47 cm	25 cm	19 cm	2,5 meter	3 cm/s (85 kg), 4 cm/s (0 kg)
Robin Mover	47 cm	25 cm	31 cm	2,5 meter	3 cm/s (85 kg), 4 cm/s (0 kg)

		
Robin	7 kg	200 kg
Robin Mover	9,5 kg	200 kg

Colour



white/grey

Warranty: 3 years

except on battery and parts of normal wear and tear (e.g. hoist tape)

Lifts per charge; Aprox 40. (user weight 85 kg, 50 cm up/down)

Voltage between battery/hoist; 24V DC

Electrical connection of charger: 230V AC 50 Hz

INVACARE® Limited
 Pencoed Technology Park - Pencoed -
 Bridgend CF35 5HZ
 Tel: +44 1 656 776 222
 Fax: +44 1656 776 220
 E-mail: uk@invacare.com
 www.invacare.co.uk
 Sales Order E-mail:
 ordersuk@invacare.com

INVACARE® Ireland Ltd
 Unit 5 - Seatown Business Campus
 Seatown Road - SWORDS - County
 Dublin
 Tel. +353 1 8107084
 Fax +353 1 8107085
 Email: Ireland@invacare.com
 www.invacare.ie

© Invacare® International Sarl 2010
 All rights reserved.

All measurement have been done under
 the Invacare® standard test conditions,
 unless quoted differently. All information
 quoted is believed to be correct at time
 of print.

Invacare® reserves the right to alter
 product specifications without prior
 consultation.

Robin - UK - 09/2010 - 1503844



Yes, you can.®